


Enhance Your Pool - One Step At A Time

Please bookmark our website, www.rehabsystems.net, as a resource for further details. Our website includes detailed descriptions, frequently asked questions and testimonials.

AquaTrek Ladder

*The AQ-5000
5-Tread Ladder*

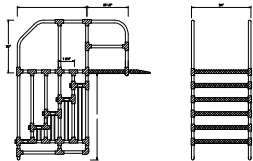


The AquaTrek Ladder is the dignified, forward-walking entry and exit to the pool that makes pool access safe and easy. Each AquaTrek Ladder is comfortable to use with wide treads and strong, metal-reinforced siderails. Stability is guaranteed due to deck mounting features, while portability is assured by hitchpins that separate the unit from the deck mounting position. Unlike other walk-in pool entry systems, the AquaTrek Ladder offers small, customizable distances between steps to guarantee a therapeutic system. To be therapeutic, a walk-in system needs to have the smallest possible distance between steps. The AquaTrek Ladder construction combines a rigid, metal skeleton with a tough, furniture-grade PVC skin. Each step tread includes a

two-tone, nonskid surface for safety. Built-in wheels on the AquaTrek Ladder allow ease of movement of the unit in and out of the pool and around the deck area. The AquaTrek Ladder is available in a standard unit rated for 600 lbs, or a heavy-duty unit to extend the standard rating. Each AquaTrek Ladder has a one-year warranty backed by a manufacturer established for over 15 years. The AquaTrek ladder is available in four, five or six step models.

AquaTrek Step


*The AQ-6000
6-Tread Step*

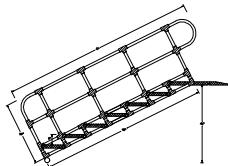

The AquaTrek Step is the next level for walk-in access to your pool. Designed as a successor unit to the AquaTrek Ladder, the AquaTrek Step is the solution for deep pool situations that demand six treads, or more, to be therapeutic. It is a small distance between steps that defines a walk-in system as therapeutic. The AquaTrek Step is available in six, seven or eight tread models.

AquaTrek Ramp

The AQ-9000 Ramp

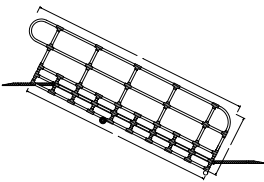


The AquaTrek Ramp is the practical solution for ramped access to a swimming pool. By using an AquaTrek Ramp system, clients receive a gentle, sloped entry where they can either walk or be wheeled into the pool. The two versions of the AquaTrek Ramp available are offered for the same value pricing. The AQ-9000 AquaTrek Ramp Standard is the ramp meant for pools less than 47-inches deep from pool floor to deck surface. The AQ-9000L AquaTrek Ramp Long is a longer length AquaTrek Ramp meant for pools where the depth from pool floor to deck surface is 48-inches or greater. The AquaTrek Ramp is designed with cross-bars to directly connect the

bottom of the two siderails, adding an extra layer of strength and durability to the unit. A 30-inch tread width on the AquaTrek Ramp assures stability, unlike wider systems where 36-inch distances between the siderails allow for more potential instability. Elevation kits are also available to help reduce the slope of an AquaTrek Ramp system in a pool. The AquaTrek Wheel Chair is the perfect companion item for use with the AquaTrek Ramp.

AquaTrek Wheel Chair

*The AQ-350-WC
Wheel Chair*

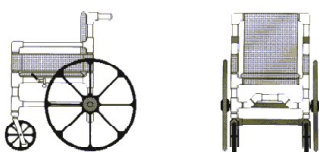

The AquaTrek Wheel Chair is totally submersible and perfect for your aquatic wheel chair needs. It is designed specifically for use in either a zero depth pool or on a ramp system like the AquaTrek Ramp. The bushings of the AquaTrek Wheel Chair are custom designed using Nylatron to prevent corrosion. A retractable footrest and soft, cushioned seat on the AquaTrek Wheel Chair provide comfort for each client. Stainless steel handrims on the rear wheels allow client mobility. A swing-arm allows for ease of mobility on and off the AquaTrek Wheel Chair. The furniture grade PVC construction of the AquaTrek Wheel Chair ensures durability for the heavy-duty needs of

the commercial aquatic environment. The options of a shoulder harness and seat belt can be purchased for the AquaTrek Wheel Chair upon request.

"The AquaTrek Ladder makes the pool much more accessible and safe for elderly members and our members with disabilities.

Thank you for such a high quality product that far outdoes our last pool ladder. "

*—Kati Tarlecky
Allentown YMCA/YWCA*

"I don't think you could destroy that (AquaTrek) in 200 years. I love this pool access system and want it! "

*—Chris Callahan
St. David's Rehab*

"It (An AquaTrek Ladder) has worked great for 13 years. "

*—Tab Ray
City of Waverly*

"It (an AQ-5600 AquaTrek Ladder) permits people to get into our pool WITHOUT using the lift chair. This is a confidence builder to people who are otherwise losing mobility. "

*—Jaimie Endresu
Mon Valley YMCA*